

# SPE - 3276

## Automated Solid Phase Extractor



EPC / PRODUCTS / APPLICATION / SOFTWARE / ACCESSORIES / CONSUMABLES / SERVICES

**Analytical Technologies Limited**

An ISO 9001 Certified Company

## ▶▶ Product Introduction

Automatic solid phase extraction series adopts highly integrated pipeline design, with multi-function sample needle, relying on precise and flexible mechanical arm to improve sample preparation efficiency and bring you convenient and easy working experience , Let your work be faster.

## ▶▶ Characteristics

- The injection needle and the extraction column are sealed by a universal plug sealing structure and equipped with H-axis pressing block to ensure good sealing.
- The suspension structure is adopted in the test tube bracket, which ensures the positioning accuracy of the manipulator to the greatest extent.
- The instrument standard router can be connected by WiFi or network line, which fully realizes the whole process of automation and unattended.
- Using modular design, only one communication line and pipeline can work between modules, and it is convenient to replace and maintain.
- Add emergency stop button, when the instrument is abnormal, manual operation can be emergency stop, to maximize user safety.
- Handwheel adjusting origin greatly saves the time of adjusting origin.

## ▶▶ Main features

**Number of channels and control mode:** The instrument host adopts modular design, supports the number of channels upgraded in the field, and supports 1/2/4/6/8 multiple channel combinations;

**Continuous processing capacity:** Maximum continuous processing capacity of 52-bit samples;

**Robot arm:** X/Y/Z three-dimensional manipulator design, X/Y axis linkage, improve operational efficiency; H axis is added to move the slider, giving consideration to flexible movement and accurate positioning;

**Sample pump:** Sampling accuracy, solvent corrosion resistance, high pressure resistance(maximum withstand voltage 1Mpa), and compatible with size and volume sampling;

**Sample needle:** 316L stainless steel, solvent resistant; unique side-bottom opening to ensure full sample loading; highly integrated function for liquid channels, gas channels, and for extraction column seals;

**Remote connection and control:** The instrument is equipped with a router that can be connected via Wi-Fi or network cable, which fully realizes the whole process of automation and unattended operation;

**Puncture sampling:** Sample needle can be selected to realize puncture sampling and avoid sample volatilization loss;

**Sample volume and receiving volume:** The whole series supports 10 - 20mL sample tube, and can support 50mL and larger specifications within two channels; the receiving volume is standard 10mL/15mL, and the two channels can be expanded by 50mL;

**Pressure sensing:** Each channel is equipped with a standard pressure sensor, push column pressure to achieve closed-loop fuzzy control, and has the functions of overpressure alarm and termination of the experiment;

**Velocity control:** The top of each extraction column is equipped with an independent sealing cap. The high pressure sealing is realized by docking with the sample needle and moving Z-axis pressure block. The solution passes through the extraction column according to the set flow rate to ensure good reproducibility and parallelism;

**Solvent management system:** Supporting 8 reagents; Supporting the hint function of insufficient residual capacity; Supporting liquid level following function;

**Infiltration and incubation:** When the active solvent or eluting solvent is loaded into the extraction column, the liquid level can be kept above the packing for a certain time to ensure full activation or elution and to ensure the ideal experimental results;

**Cleaning mode:** Supporting various cleaning modes to ensure that the inner and outer walls of sample needles are clean without cross-contamination;

**Waste discharge passage:** Support environmental protection to recover waste liquor and separate treatment of waste liquor with different properties;

**Specification of extraction column:** Compatible with conventional 1/3/6mL extraction column, other specifications of extraction column can be customized according to requirements;

**Multi-scheme extraction:** Supporting the simultaneous operation of different specifications of extraction columns and different treatment projects without waiting, improving the efficiency of continuous treatment;

**A/B column extraction:** Support multi-mode tandem column extraction, automatic identification of A/B column extraction scheme;

**Online drying/dewatering:** Supporting on-line nitrogen purge mode;

**Online Concentration:** Supporting in-situ nitrogen purge and concentration, the concentration time can be edited;

**Large volume injection:** Supports the continuous and smooth loading of more than 500 mL samples without interruption, particle wear resistance, corrosion resistance, no need for any consumables and regular maintenance; supports the continuous switching of samples, the number of samples larger than 48 bits;

**Cloud services:** Supporting upload and download programs and parameters; Supporting online technical support and exchange; Supporting remote monitoring;

**Booking starts:** After editing the experimental scheme, the start-up time can be self-booked according to the time of the analyzer;

**Count down:** After the beginning of the experiment, the countdown time is used to facilitate the control of the end time of the experiment and to match the time of the analyzer in advance;

**Intelligent self-check and alarm:** After the parameter setting is completed, the software self-checks whether there are logic errors, promptly reminds the modification; Intelligent fault alarm system, according to the fault level, chooses to continue or stop.

## ►► Technical Index

	SPE3276	SPEX00 (Customization)
Channels	4	X
Channel Independent Control	Support	
Sample pump	Precision Continuous Sample Pump	
Pump flow rate	0.1-120mL/min	
Types of solvents	8Seed, 500 mL Solvent Bottle	
Sample volume	0.1mL-20L	
Sample tube volume	10 - 20 mL,customization	
Receiving volume	10mL/15mL, other specifications customized	
Puncture sampling	Supporting	
Loading mode	X/Y/Z Three-Dimensional Manipulator	
Accuracy of manipulator	0.01mm	
Sealing technology	Docking Technology of Independent Seal Cover and Pan-plug Seal	
Extraction mode	Positive pressure extraction	
column pressure	≤80psi	
Pressure sensing	Support, Overvoltage Alarm, Threshold Editable (5-145psi)	
Column drying method	Nitrogen purge	
Waste collection	Branch collection, 3 channels	
String string	Support, A/B Pillar	
Cleaning mode	Nitrogen purge, solvent cleaning, blowing stirring cleaning and other modes are used to clean the inner and outer walls of sample needles.	
Liquid level following	Support	
Solvent Residual Capacity	Solvent deficiency hints	
Host Size (H*W*D)	680mm x 580mm x 670mm	
Voltage	220/110V, 50/60Hz	
working temperature	5-40Temperature	
Working humidity	0-90% without condensation	
power		<200W

## HPLC Servicing, Validation, Trainings and Preventive Maintenance :

**HPLC Servicing** :HPLC Servicing : We have team of service engineers who can attend to any make of HPLC promptly @the most affordable cost.

**Trainings** :We also take up preventive Maintenance to reduce downtime of HPLC's Trainings.

**AMC's/CMC** :AMC's/CMC :We offer user training both in-House and at customer sites on HPLC principles, operations, trouble-shooting.

**Validations** :Validations :We have protocols for carrying out periodic Validations as per GLP/GMP/USFDA norms.

**Instruments** :Instruments :We offer instruments/Renting Services Modules like pumps,detector etc. on Rent.



## About Analytical Technologies

Analytical Technologies is synonymous for offering technologies for doing analysis and is the Fastest Growing Global Brand having presence in at least 96 countries across the global. Analytical Technologies Limited is an ISO:9001 Certified Company engaged in Designing, Manufacturing, Marketing & providing Services for the Analytical, Chromatography, Spectroscopy, Bio Technology, Bio Medical, Clinical Diagnostics, Material Science & General Laboratory Instrumentation. Analytical Technologies, India has across the Country operations with at least 4 Regional Offices, 6 Branch Offices & Service Centers. Distributors & Channel partners worldwide.

## Our Products & Technologies




UV/VIS  
Spectro 2080+  
Double Beam



Infra FTIR



Optima Gas  
Chromatograph  
3007

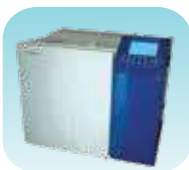

Optima Gas  
Chromatograph  
2979 Plus



Flash  
Chromatograph



Atomic Absorption  
Spectrophotometer



Liquid Particle  
Counter



Optical Emission  
Spectrophotometer



DSC/TGA



Semi Auto Bio  
Chemistry Analyzer



HEMA 2062  
Hematology  
Analyzer



Micro Plate  
Reader/Washer



URINOVA 2800  
Urine Analyzer



Total Organic  
Carbon 3800



Fully Automated  
CLIA



NOVA-2100  
Chemistry Analyzer



PCR/Gradient PCR/  
RTPCR



TOC  
Analyzer



Laser Particle  
Size Analyzer



Ion Chromatograph



Water purification  
system

## Regulatory compliances



## Corporate Social Responsibility

Analytical Foundation is a nonprofit organization (NGO) found for the purpose of:



# Analytical Foundation

1. Research & Innovation Scientist's awards/QC Professional Award : Quality life is possible by innovation only and the innovation is possible by research only, hence ANALYTICAL FOUNDATION is committed to identify such personalities for their contributions across various field of Science and Technology and awarding them yearly. To participate for award, send us your details of research / testing / publication at [Info@analyticalfoundation.org](mailto:Info@analyticalfoundation.org)

2. Improving quality of life by offering YOGA Training courses, Work shops/Seminars etc.

3. ANALYTICAL FOUNDATION aims to DETOXYFY human minds,souls and body by means of yoga, Meditation, Ayurveda, Health Care, Awards, Media, Events, Camps etc.



**HPLC Solutions   MultipleLabs   Analytical Bio-Med   Analytical Distributors   Analytical Foundation (Trust)**

**Corporate & Regd. Office:**  
Analytical House, # E67 & E68,  
Ravi Park, Vasna Road, Baroda,  
Gujarat 390 015. INDIA

T: +91 265 2253620  
+91 265 2252839  
+91 265 2252370  
F: +91 265 2254395

E: [info@hplctechnologies.com](mailto:info@hplctechnologies.com)  
[info@multiplelabs.com](mailto:info@multiplelabs.com)  
[info@analyticalgroup.net](mailto:info@analyticalgroup.net)  
[info@analyticalbiomed.com](mailto:info@analyticalbiomed.com)

W: [www.analycalgroup.net](http://www.analycalgroup.net)  
[www.hplctechnologies.com](http://www.hplctechnologies.com)  
[www.multiplelabs.com](http://www.multiplelabs.com)  
[www.ais-india.com](http://www.ais-india.com)

Sales & Support Offices:  
across the country :  
Distributors & Channel  
partners World Wide