

BGCA - 3015

Blood Gas and Chemistry Analyzer



EPCC / PRODUCTS / APPLICATION / SOFTWARE / ACCESSORIES / CONSUMABLES / SERVICES

Analytical Technologies Limited

An ISO 9001 Certified Company

www.analyticalgroup.net

Parameter	Reportable Range	Resolution	Measurement method
pH (pH units)	6.5 – 8.0	0.001	Potentiometric sensor
pO ₂ (mmHg)	10 – 700	0.1	Amperometric sensor
pCO ₂ (mmHg)	10 – 150	0.1	Potentiometric sensor
Na (mmol/L)	100 – 180	0.1	Potentiometric sensor
K (mmol/L)	2.0 – 9.0	0.01	Potentiometric sensor
Ca (mmol/L)	0.25 – 2.5	0.01	Potentiometric sensor
Cl (mmol/L)	65 – 140	0.1	Potentiometric sensor
Glu (mg/dL or mmol/L)	20 – 700/1.1 – 38.9	1/0.1	Amperometric sensor
Lac (mg/dL or mmol/L)	2.7 – 180.2/0.3 -20	0.1/0.01	Amperometric sensor
Hct (%PCV or Fraction)	10 – 75 /0.10 – 0.75	1%/0.01	Conductance sensor

►► Specification:

Throughput	Results in 1 minute after sample aspiration
Sample volume	140µl
Quality control	3 or 5 levels QC, External electronic simulator
Display	7-inch color LCD Display, 800*480
Interface	4 x USB 2.0 host, 1 x RS232, WLAN
Input device	Touch screen and barcode scanner
Power supply	100~240VAC, 50/60Hz
Battery	5000mAh rechargeable lithium-ion battery, 50 samples continuous testing
Dimensions (W*H*D)	238*153*310 mm
Weight	3.65 KG
Operation Environment	10 °C-31 °C; %RH: 25%-80%; 70 -106.6 KPa

LAN/Wi-Fi connectivity



Build-in Printer and USB Storage



Printer



Cartridge



Interface

USB Storage



Traffic Lights



Power Switch



Portable, lightweight

- Capable to run 50 samples with fully-charged battery
- Diagnose at the point of care, patient side, out in the field or exam room

Easy, Quick and Convenient

- Zero maintenance
- Minimize hands on time and training requirement
- Auto-sampling
- Report is ready within one minute after sample aspiration

Accurate and Reliable

- Innovative microchip liquid control technology and micro-sensor multifunction membrane technology
- High sensitivity and accuracy
- Calibrator, Control and Electronic simulator to TRIPLE guarantee the RESULTS!

Flexible Data Management

- Acquire patients' information via barcode/QR code
- Up to 10,000 patient data storage
- 4 USB ports and LAN/WiFi for data management
- Optional data management software

Multi-parameter cartridge

- Multi popular time-sensitive parameters come in one cartridge, including ABG, electrolytes and metabolites
- Single-use cartridge avoids contamination
- Room temperature storage with long shelf-life



Portable and lightweight

Touch screen with built-in multimedia tutorials



Variety of test Cartridges

	pH	pCO ₂	pO ₂	Na	K	Cl	Ca	Hct	Glu	Lac
BG3	=====	=====	=====							
BG8	=====	=====	=====	=====	=====	=====	=====	=====		
BG4	=====	=====	=====							=====
BG9	=====	=====	=====	=====	=====	=====	=====	=====	=====	
BG10	=====	=====	=====	=====	=====	=====	=====	=====	=====	=====
BC4				=====	=====	=====	=====	=====		

Calculated values: ch+, HCO₃-act, HCO₃-std, BE(ecf), BE(B), BB(B), ctCO₂, sO₂(est), Ca⁺⁺(7.4), AnGap, tHb(est), pO₂(A-a), pO₂(a/A), Ri, pO₂/FIO₂, ch+(T),pH(T), pCO₂(T), pO₂(A-a)(T), pO₂(a/A)(T), RI(T), pO₂(T)/FIO₂, mOsm

HPLC Servicing, Validation, Trainings and Preventive Maintenance :

HPLC Servicing :HPLC Servicing : We have team of service engineers who can attend to any make of HPLC promptly @the most affordable cost.

Trainings :We also take up preventive Maintenance to reduce downtime of HPLC's Trainings.

AMC's/CMC :AMC's/CMC :We offer user training both in-House and at customer sites on HPLC principles, operations, trouble-shooting.

Validations :Validations :We have protocols for carrying out periodic Validations as per GLP/GMP/USFDA norms.

Instruments :Instruments :We offer instruments/Renting Services Modules like pumps,detector etc. on Rent.



About Analytical Technologies

Analytical Technologies is synonymous for offering technologies for doing analysis and is the Fastest Growing Global Brand having presence in at least 96 countries across the global. Analytical Technologies Limited is an ISO:9001 Certified Company engaged in Designing, Manufacturing, Marketing & providing Services for the Analytical, Chromatography, Spectroscopy, Bio Technology, Bio Medical, Clinical Diagnostics, Material Science & General Laboratory Instrumentation. Analytical Technologies, India has across the Country operations with at least 4 Regional Offices, 6 Branch Offices & Service Centers. Distributors & Channel partners worldwide.

Our Products & Technologies



UV/VIS
Spectro 2080+
Double Beam



Infra FTIR



Optima Gas
Chromatograph
3007



Optima Gas
Chromatograph
2979 Plus



Flash
Chromatograph




Atomic Absorption
Spectrophotometer



Liquid Partical
Counter



Optical Emission
Spectrophotometer

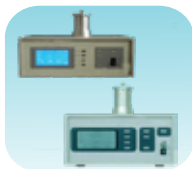

DSC/TGA



Semi Auto Bio
Chemistry Analyzer



HEMA 2062
Hematology
Analyzer



Micro Plate
Reader/Washer



URINOVA 2800
Urine Analyzer



Total Organic
Carbon 3800



Fully Automated
CLIA



NOVA-2100
Chemistry Analyzer



PCR/Gradient PCR/
RTPCR



TOC
Analyzer



Laser Particle
Size Analyzer



Ion Chromatograph



Water purification
system

Regulatory compliances



Corporate Social Responsibility

Analytical Foundation is a nonprofit organization (NGO) found for the purpose of:



Analytical Foundation

1. Research & Innovation Scientist's awards/QC Professional Award : Quality life is possible by innovation only and the innovation is possible by research only, hence ANALYTICAL FOUNDATION is committed to identify such personalities for their contributions across various field of Science and Technology and awarding them yearly. To participate for award, send us your details of research / testing / publication at Info@analyticalfoundation.org

2. Improving quality of life by offering YOGA Training courses, Work shops/Seminars etc.

3. ANALYTICAL FOUNDATION aims to DETOXYFY human minds,souls and body by means of yoga, Meditation, Ayurveda, Health Care, Awards, Media, Events, Camps etc.

Reach us @



HPLC Solutions MultipleLabs Analytical Bio-Med Analytical Distributors Analytical Foundation (Trust)

Corporate & Regd. Office:
Analytical House, # E67 & E68,
Ravi Park, Vasna Road, Baroda,
Gujarat 390 015. INDIA

T: +91 265 2253620
+91 265 2252839
+91 265 2252370
F: +91 265 2254395

E: info@hplctechnologies.com
info@multiplelabs.com
info@analyticalgroup.net
info@analyticalbiomed.com

W: www.analycalgroup.net
www.hplctechnologies.com
www.multiplelabs.com
www.ais-india.com

Sales & Support Offices:
across the country :
Distributors & Channel
partners World Wide

Note : Company reserves rights to add/delete/modify the contents / technical specifications of the catalogue without prior notice.