



Cell Saver CS-3025

Autologous Blood Recovery System

Performance you can count on –
because patients are counting on you.



EPC / PRODUCTS / APPLICATION / SOFTWARE / ACCESSORIES / CONSUMABLES / SERVICES

Analytical Technologies Limited

An ISO 9001 Certified Company

www.analyticalgroup.net

Cell Saver

You know you're at your best when everything is working as it should


When the surgeon is focused on the surgery. And the patient is getting what's best—their own high-quality blood. Your standards are high. That's why you demand excellence.

We've evolved cell salvage to a higher level...yours.

The Cell Saver system provides hospitals with an easy-to-use, reliable way to recover and deliver back a patient's own high-quality blood during medium to high blood loss procedures, including cardiac, orthopedic, trauma, transplant, vascular, and OBGYN surgeries. The complete system is designed to be intuitive and flexible for optimal performance and to support exceptional patient care.


What a higher level of performance looks like:

- ✓ Everything running smoothly.
- ✓ A focused surgeon.
- ✓ And a well cared for patient.
- ✓ All supported by a state-of-the-art system.

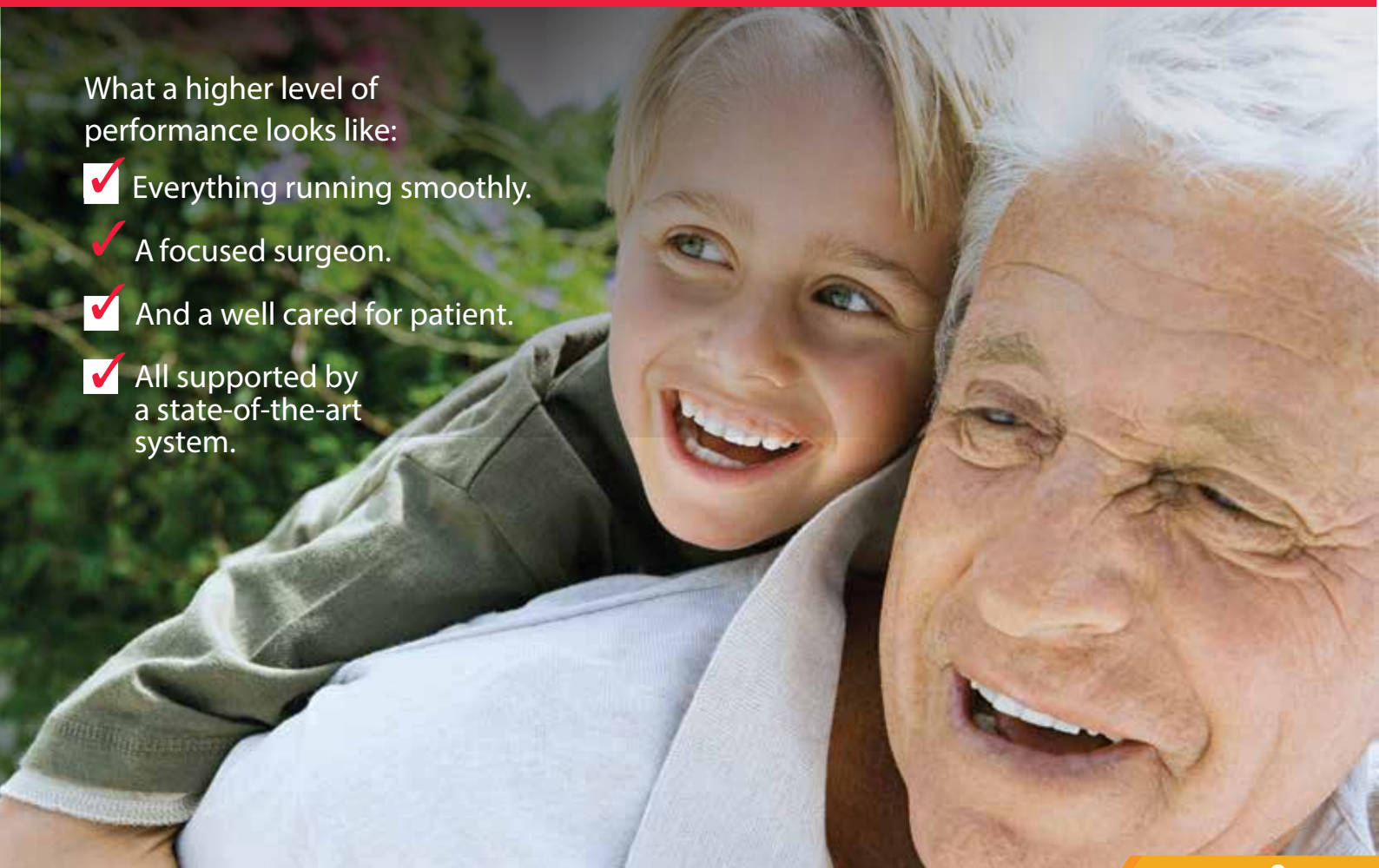

Focus on quality and ease of use

The Cell Saver provides patients with their own fresh, washed blood, preventing an allogeneic transfusion while helping to reduce costs. Its integrated SmartSuction is gentle on red blood cells (RBC's), preserving them for transfusion. And, its fat reduction protocol is designed for optimal performance, while still offering custom settings that allow for modifications during a procedure.

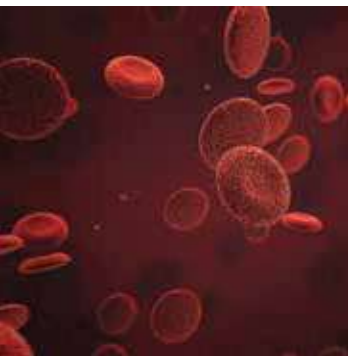

This results in:


- High hematocrit with excellent RBC recovery—providing more red cells for transfusion.
- Optimal washout of undesirable components.
- Exceptional fat removal performance.
- Maximized patient safety—with autotransfusion and sequestration.

High Quality Red Blood Cells (RBC's) ¹

• Hematocrit (Hct) level	~ 60%
• RBC recovery	~ 94%
• Washout performance	
✓ Albumin removal	~ 98%
✓ Free Hemoglobin	~ 99%
✓ Heparin removal	~ 99%
✓ Potassium	~ 96%
• Fat Removal > 99% ²	

Integrated SmartSuction

- Protects RBC's while maintaining optimal suction.
- Significantly reduce hemolysis compared to constant pressure. ³
- More accurate and reliable than wall suction.

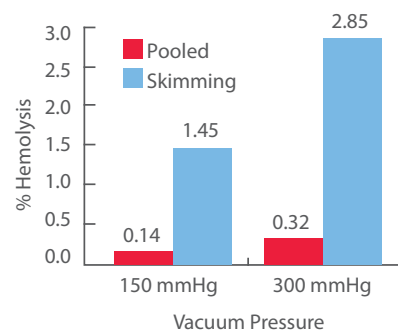

Gregoretti et al. showed that greater hemolysis is associated with higher suction forces. Hemolysis is even greater when air mixes with blood. ⁴

Cell Saver

Fast, reliable, and intuitive


The system is fast, reliable, and intuitive, incorporating state-of-the-art software enhancements.

You'll have the flexibility to define and save password-protected protocols to meet unique surgical needs. The system's touch screen instructions and messaging are easy to read and navigate, complete with on-screen help. Its status beacon is highly visible, from any direction, and the barcode scanner captures disposable, patient, and surgical procedure information.



User friendly, full color touch screen

Intuitive and simple programming

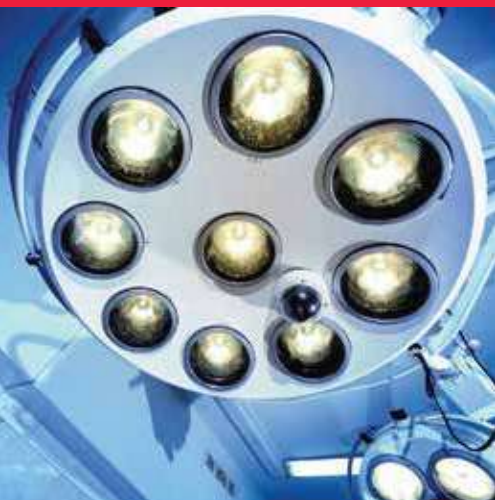

- Easy to navigate touch screen with color-coded messaging and corresponding status beacon that can be seen throughout the OR.
- Fast start-up time with suction options available in seconds.
- On-screen help is always available.
- Minimal training necessary to keep staff proficient—meaning more time focused on the patient.

- One touch operation with flexibility to customize and save settings for your unique needs.
- Ability to password protect settings groups.
- Custom alarm tones and volume controls.
- Preset procedure record information for ease of entry.
- Flexibility to define settings group for system start-up.



Easy to transport between ORs

- Smaller footprint, lightweight, and portable.
- Multiple storage options.
- Durable casters with easy
- Components fold away for easy transport.

Fast, easy, disposable set-up

- Cell salvage and sequestration sets are color-coded for ease of use.
- Single barcode scanning for all disposable sets and ease of data entry.
- Multiple bowl volume options and packaging configurations.



70 mL



125 mL



225 mL

Cell Saver

Secure, data collection, management and transfer capabilities

- Minimize human errors and associated costs through accurate data collection and transfer.
- Capture and store up to 100 procedures.
- Track and trend procedure data effortlessly with procedure records downloaded directly to your PC.
- Technical data downloads can be accessed for analysis and best practice sharing.



Meeting your product performance needs of today, while planning for tomorrow.

Analytical meets your clinical performance needs of today, while planning for tomorrow. The Cell Saver software architecture and system capabilities allow for continued advancement—all while providing secure data collection and transfer.

This results in:

- The best care for patients.
- A reduced need for allogeneic transfusions.
- A commitment to high-quality product engineering, design, manufacturing, and reliability.
- A platform that allows for quick software upgrades.



Ordering Information

Description	List Number	Qty per Case
Cell Saver System	CS3025	1
Cell Saver Cart	CS-EA-1000	1
Cell Saver Fast Pack (225 mL, 205, 208, HAR-A-1000)	CS-FP-225V	4
Cell Saver Fast Pack (225 mL, 220, 208, HAR-A-1000)	CS-FP-225F	4
Cell Saver Fast Pack (125 mL, 205, 208, HAR-A-1000)	CS-FP-125V	4
Cell Saver Processing Set (225 mL)	CS-P-225	8
Cell Saver Processing Set (125 mL)	CS-P-125	8
Cell Saver Processing Set (70 mL)	CS-P-70	8
Sequestration Kit	CS-SQ-1000	8
Collection reservoir w/150 µ filter, 3 L	00205-00	4
Collection reservoir w/20 µ filter, 3 L	00220-00	4
Cell Saver Aspiration & Anticoagulation Line	00208-00	20
SmartSuction Aspiration & Anticoagulation Line	ATL-A-1003	10
SmartSuction Filtered Vacuum Line, non-sterile	ATL-A-1000	10
RBC Bag, 1000 mL	00245-00	40
Waste Bag, 10 L	CS-B-1000	10
70 mL Bowl Chuck Adapter	107581-00	1
40 µm Microaggregate Blood Transfusion Filter	SQ40S	40
40 µm Microaggregate Blood Transfusion Filter with attached straight administration set	SQ40SK	16
40 µm Microaggregate Blood Transfusion Filter with attached Y-type administration set	SQ40SY	16
LipiGuard Filter for Salvaged Blood Stand Alone (no administration set attached)	SB1E	20
RC Leukocyte Reduction Filter for Intraoperately Salvaged Washed Blood Stand Alone (no administration set attached)	RS1	20
Post-op Line	01500-00	20
Extended Post-op Line	01550-00	20

Specifications

Device only	Height: 16.5 in (41.9 cm); Width: 11.75 in (29.8 cm); Depth: 21.5 in (54.6 cm)
Device with cart	Height: 72 in (182.9 cm) poles extended; 48 in (121.9 cm) poles down; Width: 21 in (53.3 cm); Depth: 26.5 in (67.3 cm)
Weight of device	< 56 lbs (25 kg)
Weight of cart	< 40 lbs (18 kg)
Voltage	100–240 VAC ±10%
Pump speed	0–1000 mL/min (adjustable)
Centrifuge speed	Latham bowl: 5650 rpm; Blow-molded bowl: 7500 rpm; Sequestration: 2050-5650 rpm (adjustable)
Fuse rating	100–240V: T3.15A250V
Operating frequency	47–63 Hz
Laser classification	Class 3R laser product
Power cord length	16 ft (4.9 m), 2 wire ground (earth) connection, 3 prong hospital grade

Analytical Partner In Practice

We offer access to procedural data and troubleshooting to support best practices, dedicated field services—including clinical specialist teams, knowledgeable sales teams, and award-winning customer service.

This results in:

- A committed partner—starting with implementation and training and supporting all ongoing needs.
- Ensured staff proficiency during full initial implementation.
- A hotline for customer support.



▶▶ **HPLC Servicing, Validation, Trainings and Preventive Maintenance :**

- HPLC Servicing :HPLC Servicing : We have team of service engineers who can attend to any make of HPLC promptly @the most affordable cost.
- Trainings :We also take up preventive Maintenance to reduce downtime of HPLC's Trainings.
- AMC's/CMC :AMC's/CMC :We offer user training both in-House and at customer sites on PLC principles, operations, troubleshooting.
- Validations :Validations :We have protocols for carrying out periodic Validations as per GLP/ GMP/USFDA norms.
- Instruments :Instruments :We offer instruments / Renting Services Modules like pumps, detector etc. on Rent.



▶▶ **About Analytical Technologies**

Analytical Technologies is synonymous for offering technologies for doing analysis and is the Fastest Growing Global Brand having presence in at least 96 countries across the globe. Analytical Technologies Limited is an ISO :9001 Certified Company engaged in Designing, Manufacturing, Marketing & providing Services for the Analytical, Chromatography, Spectroscopy, Bio Technology, Bio Medical, Clinical Diagnostics, Material Science & General Laboratory Instrumentation. Analytical Technologies, India has across the Country operations with at least 4 Regional Offices, 6 Branch Offices & Service Centers. Distributors & Channel partners worldwide.

▶▶ **Our Products & Technologies**

						
UV/VIS Spectro 2060+ Double Beam	FTIR	Gas Chromatograph 3000	Gas Chromatograph 2979 Plus	Flash Chromatograph	Atomic Absorption Spectrophotometer	Liquid Partical Counter
						
Optical Emission Spectrophotometer	DSC/TGA	NOVA 2020 plus Automated Bio Chemistry Analyzer	HEMA 2020 Hematology Analyzer	Micro Plate Reader/Washer	URINOVA 2800 Urine Analyzer	Total Organic Carbon
						
Fully Automated CLIA	NOVA Basic Semi-Auto Chemistry Analyzer	PCR/Gradient PCR/ RTPCR	Blood Gas Analyzer	Random access Analyzer for immunoassay Proteins & clinical chemistry	Semen Analyzer	Water purification system

▶▶ **Regulatory compliances**



▶▶ **Corporate Social Responsibility**

Analytical Foundation is a Nonprofit Organization (NGO) found for the purpose of:



Analytical
Foundation

1. Research & Innovation Scientist's awards / QC Professional Award : Quality life is possible by innovation only and the innovation is possible by research only, hence ANALYTICAL FOUNDATION is committed to identify such personalities for their contributions across various field of Science and Technology and awarding them yearly. To participate for award, send us your details of research / testing / publication at info@analyticalfoundation.org

2. Improving quality of life by offering YOGA Training courses, Work shops / Seminars etc.

3. ANALYTICAL FOUNDATION aims to DETOXYFY human minds, souls and body by means of Yoga, Meditation, Ayurveda, Health Care, Awards, Media, Events, Camps etc.

▶▶ **Reach us @**



HPLC Solutions	MultipleLabs	Analytical Bio-Med	Analytical Distributors	Analytical Foundation (Trust)
Corporate & Regd. Office: Analytical House, # E67 & E68, Ravi Park, Vasna Road, Baroda, Gujarat 390 015. INDIA	T: +91 265 2253620 +91 265 2252839 +91 265 2252370 F: +91 265 2254395	E: info@hplctechnologies.com info@multiplelabs.com info@analyticalgroup.net info@analyticalbiomed.com	W. www.ais-india.com www.analyticalgroup.net www.hplctechnologies.com www.multiplelabs.com	Sales & Support Offices: across the country : Distributors & Channel partners World Wide

Note : Company reserves rights to add/delete/modify the contents / technical specifications of the catalogue without prior notice.