

OxygenC - 3801

Home/Transport Oxygen Concentrator



EPCC / PRODUCTS / APPLICATION / SOFTWARE / ACCESSORIES / CONSUMABLES / SERVICES

Analytical Technologies Limited

An ISO 9001 Certified Company


www.analyticalgroup.net

►► Features

1. PSA Technology
 2. Compact, light weight and easy to carry
 3. Perfect structure and Portable Design
 4. Energy efficient and low power consumption (25W)
 5. Can be used in car voltage 12V.
 6. HEPA filter and bacteria filter can filter most of impurity, bacteria and PM2.5
- Atomizing and Timing function

►► Specification

Name:	Home/Transport Oxygen Concentrator
Model:	OxygenC - 3801
Oxygen making mode:	VPESA mode
Oxygen concentration:	45% ±3%
Oxygen Flow rate:	1000ML/min
Inhalation way:	Nasal
Spray rate:	≥0.2ML/min
Working air flow:	≥6.7LPM
Working pressure range:	≥83kPa
Power:	30W
Negative ions:	300,000/ cm ³ (electron reflection)
Weight:	2.65kg
Voltage:	DC 12V
Noise:	≤60dB 1 meter away
Dimension:	15cm x 24cm x 42cm

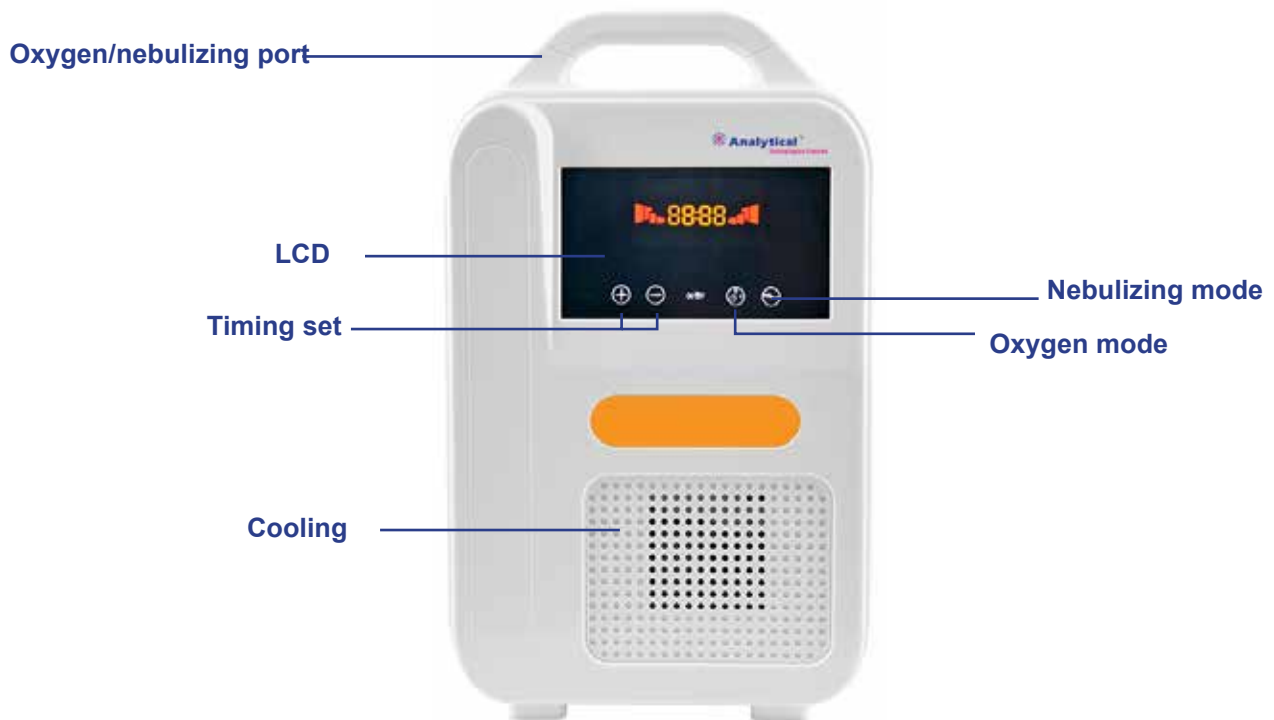

►► Accessories



Nasal Breathing Tube



Nebulizing Accessories



Adapter



Vehicle Power Cord

►► Usage

1. Plug the power adapter into the power outlet and touch the “ON/OFF” to turn on the machine.
2. Oxygen absorption: put on the nasal oxygen tube and touch the oxygen button.
Turn the knob on the top to change flow rate.
3. Atomizing mode: pour the medicine into the atomizing can and put on the atomizer,
touch the atomizing button to start.
4. Timing: touch “+” and “-” to set a fixed time.
After use, turn off the power switch and unplug the adapter.

►► Maintenance Method

1. Host maintenance method

The dirt on the main unit can be wiped gently with a soft rag.

2. Nasal oxygen tube maintenance method

Clean the nasal oxygen tube with water. Gently wipe the oxygen outlet and the nasal oxygen tube with an alcohol-based paper towel.

The machine has intermittent exhaust sound when use,
the interval is about 3second

HPLC Servicing, Validation, Trainings and Preventive Maintenance :

HPLC Servicing :HPLC Servicing : We have team of service engineers who can attend to any make of HPLC promptly @the most affordable cost.

Trainings :We also take up preventive Maintenance to reduce downtime of HPLC's Trainings.

AMC's/CMC :AMC's/CMC :We offer user training both in-House and at customer sites on HPLC principles, operations, trouble-shooting.

Validations :Validations :We have protocols for carrying out periodic Validations as per GLP/GMP/USFDA norms.


Instruments :Instruments :We offer instruments/Renting Services Modules like pumps,detector etc. on Rent.



About Analytical Technologies

Analytical Technologies is synonymous for offering technologies for doing analysis and is the Fastest Growing Global Brand having presence in at least 96 countries across the global. Analytical Technologies Limited is an ISO:9001 Certified Company engaged in Designing, Manufacturing, Marketing & providing Services for the Analytical, Chromatography, Spectroscopy, Bio Technology, Bio Medical, Clinical Diagnostics, Material Science & General Laboratory Instrumentation. Analytical Technologies, India has across the Country operations with at least 4 Regional Offices, 6 Branch Offices & Service Centers. Distributors & Channel partners worldwide.

Our Products & Technologies

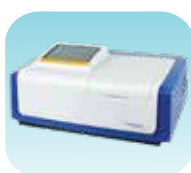

UV/VIS
Spectro 2080+
Double Beam



Infra FTIR



Optima Gas
Chromatograph
3007



Optima Gas
Chromatograph
2979 Plus



Flash
Chromatograph



Atomic Absorption
Spectrophotometer



Liquid Partical
Counter



Optical Emission
Spectrophotometer



DSC/TGA



Semi Auto Bio
Chemistry Analyzer



HEMA 2062
Hematology
Analyzer



Micro Plate
Reader/Washer



URINOVA 2800
Urine Analyzer



Total Organic
Carbon 3800



Fully Automated
CLIA



NOVA-2100
Chemistry Analyzer



PCR/Gradient PCR/
RTPCR



TOC
Analyzer



Laser Particle
Size Analyzer



Ion Chromatograph



Water purification
system

Regulatory compliances



Corporate Social Responsibility

Analytical Foundation is a nonprofit organization (NGO) found for the purpose of:



Analytical Foundation

1. Research & Innovation Scientist's awards/QC Professional Award : Quality life is possible by innovation only and the innovation is possible by research only, hence ANALYTICAL FOUNDATION is committed to identify such personalities for their contributions across various field of Science and Technology and awarding them yearly. To participate for award, send us your details of research / testing / publication at Info@analyticalfoundation.org

2. Improving quality of life by offering YOGA Training courses, Work shops/Seminars etc.

3. ANALYTICAL FOUNDATION aims to DETOXYFY human minds,souls and body by means of yoga, Meditation, Ayurveda, Health Care, Awards, Media, Events, Camps etc.

Reach us @



 **Analytical**®
Technologies Limited

HPLC Solutions MultipleLabs Analytical Bio-Med Analytical Distributors Analytical Foundation (Trust)

Corporate & Regd. Office:
Analytical House, # E67 & E68,
Ravi Park, Vasna Road, Baroda,
Gujarat 390 015. INDIA

T: +91 265 2253620
+91 265 2252839
+91 265 2252370
F: +91 265 2254395

E: info@hplctechnologies.com
info@multiplelabs.com
info@analyticalgroup.net
info@analyticalbiomed.com

W. www.analyticalgroup.net
www.hplctechnologies.com
www.multiplelabs.com
www.ais-india.com

Sales & Support Offices:
across the country :
Distributors & Channel
partners World Wide