

# Walk-In-Humidity Chamber



## Technical Specification

**Temperature Range:** 20°C - 60 °C

**Humidity Range:** 40%RH - 95%RH

**Accuracy:** +0.2 °C / +2.0 % RH

**Uniformity:** +1.0 °C / +3.0 % RH

**Test suitable for:** 25°C & 60%RH, 30°C & 65%RH,  
40°C & 75%RH, 30°C & 75%RH

Specially designed chambers for Low Humidity conditions, such as 25°C & 40%RH, 40°C & 25%RH, 30°C & 35%RH.

## Features

### Construction

Chamber interior is fabricated of stainless steel with dull finish & exterior is of stainless steel / G.I. Powder coated. The inner working chamber is water & vapour tight high-grade stainless steel. All seams are arc-welded & passivated for strength & durability. Modular type prefabricated metal sectional panels are lightweight for easy handling and also fit through standard doorways, easily transported and assembled /disassembled for relocating the room.

### Full View Glass Door

Full door is of toughened glass with polysulphide to view the samples without disturbing the test conditions & human safety.

### Trays & Racks

Heavy duty solid stainless steel duly electro polished wiremesh type trays are provided which are easy to clean, time & effort saving, adjustable for convenience, along with detachable racks with adjustable angles.

### Scanner

For Validation purpose Temp / Humidity multipoint Scanner can be provided with printer & PC interface facility.

### 21 CFR PART-11 Compliance

Log data, Event data, Audit Trail data, E-Records & E-signatures, Graphical analysis & Data acquisition.

### Stand by Refrigeration & Humidity System

Two steam generators and two water reserve tanks are provided. In case of any system failure or breakdown the stand by system is switched over automatically without disturbing the test condition.

## Optional Features

GSM TECHNOLOGY

HOOTER SYSTEM

EXTRA TRAY



Stand by Refrigeration system



Rack



Control system



## Walk-In-Humidity Chamber

### Sizes:

Models	AT*-WH600	AT*-WH800	AT*-WH1000	AT*-WH1200	AT*-WH1600	AT*-WH1800
Features	6000 Ltrs	8000 Ltrs	10000 Ltrs	12000 Ltrs	16000 Ltrs	18000 Ltrs

#### Interior Dimension

	AT*-WH600	AT*-WH800	AT*-WH1000	AT*-WH1200	AT*-WH1600	AT*-WH1800
Width in mm	2000	2000	2000	2000	2000	4000
Depth in mm	1500	2000	2500	3000	4000	3000
Height in mm	2000	2000	2000	2000	2000	2000

#### Exterior Dimension

	AT*-WH600	AT*-WH800	AT*-WH1000	AT*-WH1200	AT*-WH1600	AT*-WH1800
Width in mm	2160	2160	2160	2160	2160	4160
Depth in mm	2010	2510	3010	3510	4510	3510
Height in mm	2190	2190	2190	2190	2190	2190

#### Racks and Trays

	AT*-WH600	AT*-WH800	AT*-WH1000	AT*-WH1200	AT*-WH1600	AT*-WH1800
No of racks	4	4	6	6	8	10
No of Trays	20	20	30	30	40	50
Trays Size in mm	600x600x30	600x900x30	700x600x30	900x600x30	900x600x30	900x600x30
Type	Wiremesh Type					

#### Power Supply

	AT*-WH600	AT*-WH800	AT*-WH1000	AT*-WH1200	AT*-WH1600	AT*-WH1800
Power Supply	Three Phase, 440V 50Hz					
Power Rating	7KW	8KW	8KW	8KW	9KW	15KW
Warranty	12 Months					

Other Sizes are also available as on requirement

Models	Material of Construction	
	Enterior	Exterior
AT1	S.S.304 Dull Finish	G.I. Powder Coated
AT2	S.S.304 Dull Finish	S.S.304 Dull Finish
AT3	S.S.316 Dull Finish	S.S.304 Dull Finish

#### NOTE:

AT\* :- AT1 / AT2 / AT3

WALK-IN BOD INCUBATOR & WALK-IN COLD CHAMBER are also available on requirement.

▶▶ We Provide EPCC(Engineering, Procurement, Construction and Consultancy) Solutions For Different Areas.

such as,

**Food Safety**

**Clinical & Medical Diagnostics**

**Environment & Pollution Control**

**Genetic**


**Mobile Lab**

**Pharma Biotech**

**Soil testing Lab**

**Water Testing**

**Forensic science**

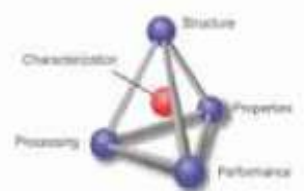

## About Analytical Technologies

Analytical Technologies is synonymous for offering technologies for doing analysis and is the Fastest Growing Global Brand having presence in at least 96 countries across the globe. Analytical Technologies Limited is an ISO :9001 Certified Company engaged in Designing, Manufacturing, Marketing & providing Services for the Analytical, Chromatography, Spectroscopy, Bio Technology, Bio Medical, Clinical Diagnostics, Material Science & General Laboratory Instrumentation. Analytical Technologies, India has across the Country operations with at least 4 Regional Offices, 6 Branch Offices & Service Centers. Distributors & Channel partners worldwide.



▶▶ Regulatory compliances



▶▶ Corporate Social Responsibility

Analytical Foundation is a Nonprofit Organization (NGO) found for the purpose of:



TM

**Analytical**  
**Foundation**

1. Research & Innovation Scientist's awards / QC Professional Award : Quality life is possible by innovation only and the innovation is possible by research only, hence ANALYTICAL FOUNDATION is committed to identify such personalities for their contributions across various field of Science and Technology and awarding them yearly. To participate for award, send us your details of research / testing / publication at [info@analyticalfoundation.org](mailto:info@analyticalfoundation.org)

2. Improving quality of life by offering YOGA Training courses, Work shops / Seminars etc.

3. ANALYTICAL FOUNDATION aims to DETOXYFY human minds, souls and body by means of Yoga, Meditation, Ayurveda, Health Care, Awards, Media, Events, Camps etc.

▶▶ Reach us @



HPLC Solutions   MultipleLabs   Analytical Bio-Med   Analytical Distributors   Analytical Foundation (Trust)

**Corporate & Regd. Office:**  
Analytical House, # E67 & E68,  
Ravi Park, Vasna Road, Baroda,  
Gujarat 390 015. INDIA

T: +91 265 2253620  
+91 265 2252839  
+91 265 2252370  
F: +91 265 2254395

E: [info@hplctechnologies.com](mailto:info@hplctechnologies.com)  
[info@multiplelabs.com](mailto:info@multiplelabs.com)  
[info@analyticalgroup.net](mailto:info@analyticalgroup.net)  
[info@analyticalbiomed.com](mailto:info@analyticalbiomed.com)

W. [www.ais-india.com](http://www.ais-india.com)  
[www.analycalgroup.net](http://www.analycalgroup.net)  
[www.hplctechnologies.com](http://www.hplctechnologies.com)  
[www.multiplelabs.com](http://www.multiplelabs.com)

Sales & Support Offices:  
across the country :  
Distributors & Channel  
partners World Wide

Note : Company reserves rights to add/delete/modify the contents / technical specifications of the catalogue without prior notice